


O-arm™

Surgical Imaging System

VERSION 4.2
FEATURE UPDATE

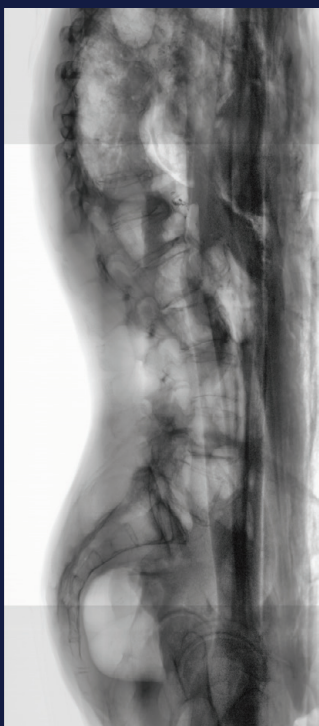

ASSESS > AFFIRM > ALIGN

O-arm™ O2 system continues to provide clinicians with advanced intra-operative tools for insightful surgical decision making through imaging.

Version 4.2 offers powerful imaging options, improved IT security functionality and features to improve the overall O-arm™ surgical imaging workflow.

Medtronic

V.4.2 FEATURES

2D LONG FILM

Measure, verify and reconcile

your patient's spinal anatomy with unique 2D long-film capability. With the combination of the robotic gantry, flat panel imaging, specialized collimator and software algorithms, you can capture the spine in both antero-posterior and lateral long films.

USER ACCESS ENHANCEMENTS

Utilize additional tools to enhance your IT security measures including:



- User sign-in credentialing using Lightweight Directory Access Protocol (LDAP)
- Audit capabilities related to protected health information (PHI), user management and system configuration
- Encryption of PHI on O-arm™ viewing station

GANTRY & ROTOR TILT DISPLAY

The gantry tilt and gantry rotor angle display enables O-arm™ imaging system users **to view the angle of the gantry** and the relative angle/location of the X-ray source & detector on the O-arm™ pendant.



Image artifacts and gray scale intensity may vary across the 2D Long Film images particularly near the blend zones. A linear artifact through the center of the 2D long film image may occur which is a result of the gantry skin seam and is normal to the design. As with other X-ray based image devices, consider best practices to minimize imaging anatomy off the plane-of-focus.



Medtronic

826 Coal Creek Circle
Louisville, CO 80027 USA
Toll Free: 877-242-9504
Telephone: 720-890-3200
Fax: 720-890-3500

CAUTION: Federal Law (USA) restricts these devices for sale by or on the order of a physician. Refer to product instruction manual/package insert for instructions, warning, precautions and contraindications.

For information on Indications, Safety, and Warnings, call Medtronic at **(877) 242-9504**, or visit Medtronic's website at www.medtronicneurosurgery.com.

medtronic.com

UC202103951 EN © 2020 Medtronic. All rights reserved. Medtronic, Medtronic logo and Further, Together are trademarks of Medtronic. All other brands are trademarks of a Medtronic company.